


Configuration FORTESSA X20

Blue Laser (488 nm)

↻ 477 **PerCP** 678 ↻
↻ 477 **PerCP-Cy5.5** (↻ 675) 696 ↻
↻ 477 **PerCP-Vio700** 704 ↻
↻ 477 **PerCP-eFluo710** 710 ↻
↻ 488 **BB700** 700 ↻

SSC

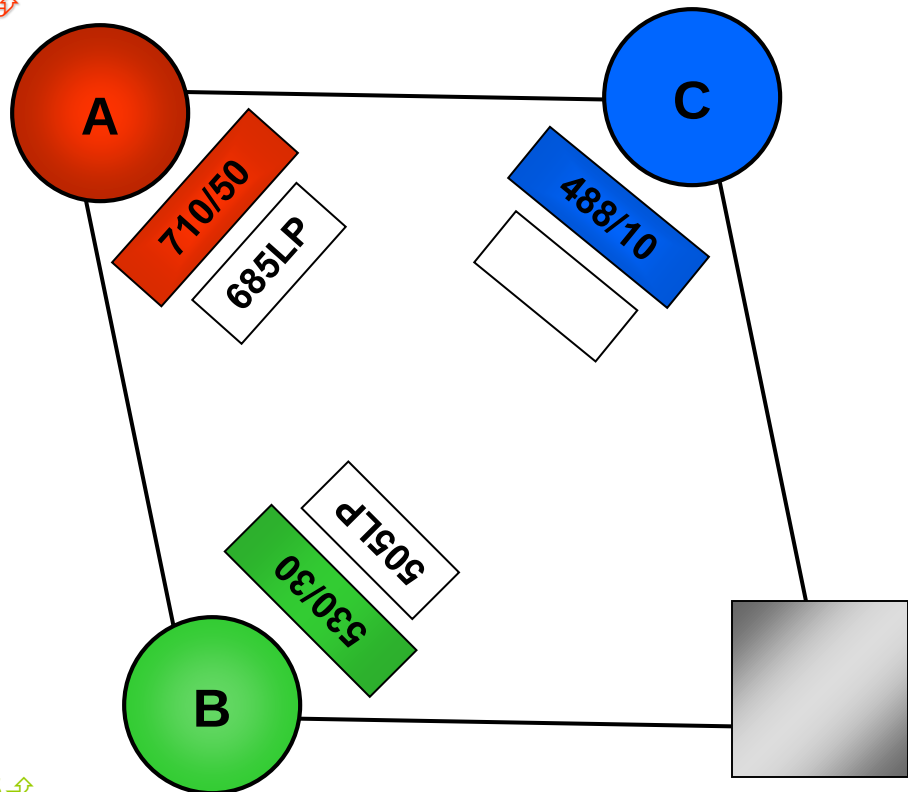


↻ 489 **GFP** 510 ↻
↻ 490 **BB515** 515 ↻
↻ 500 **Lysotracker Green** 510 ↻
↻ 493 **Calcein** 514 ↻
↻ 495 **FITC** 519 ↻
↻ 499 **Alexa Fluor 488** 520 ↻
↻ 495 **CFSE** 525 ↻
↻ 496 **VioBright FITC** 522 ↻
↻ 488 **VioBright 515** 515 ↻
↻ 488 **Vio 515** 515 ↻
↻ 490 **BD Horizon BB515** 515 ↻
↻ 490 **Zombie Green** 515 ↻
↻ 498 **BD FVS 520** 521 ↻
↻ 488 **Kiravia Blue 520** 520 ↻

Configuration FORTESSA X20

UV Laser (355 nm)

↗ 348 BUV737 737 ↗

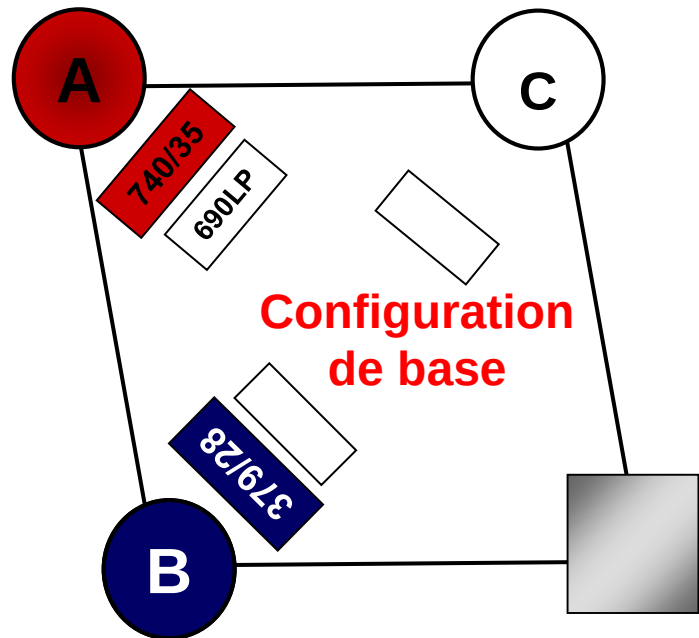


↗ 348 BUV395 395 ↗

Calcium

↗ 345 Indo 1 libre 485 ↗

↗ 382 Zombie Aqua FV 515 ↗

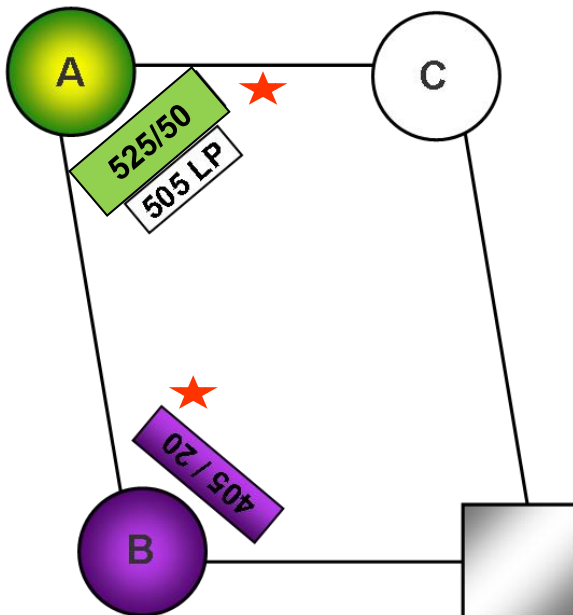


↗ 325 Indo 1 lié 400 ↗

Configuration FORTESSA X20

Yellow Green Laser (561 nm)

↻ 565 **PE-Vio615** 519 ↻
↻ 565 **PE-eFluor610** 607 ↻
↻ 566 **PE-Dazzle594** 610 ↻
↻ 496 **PE-CF594** 612 ↻
↻ 566 **PE-Texas-Red (ECD)**(↻ 595) 616 ↻
↻ 566 **PE-Alexa Fluor 610** (↻ 610) 627 ↻

↻ 566 **PE-Cy5.5** (↻ 675) 693 ↻

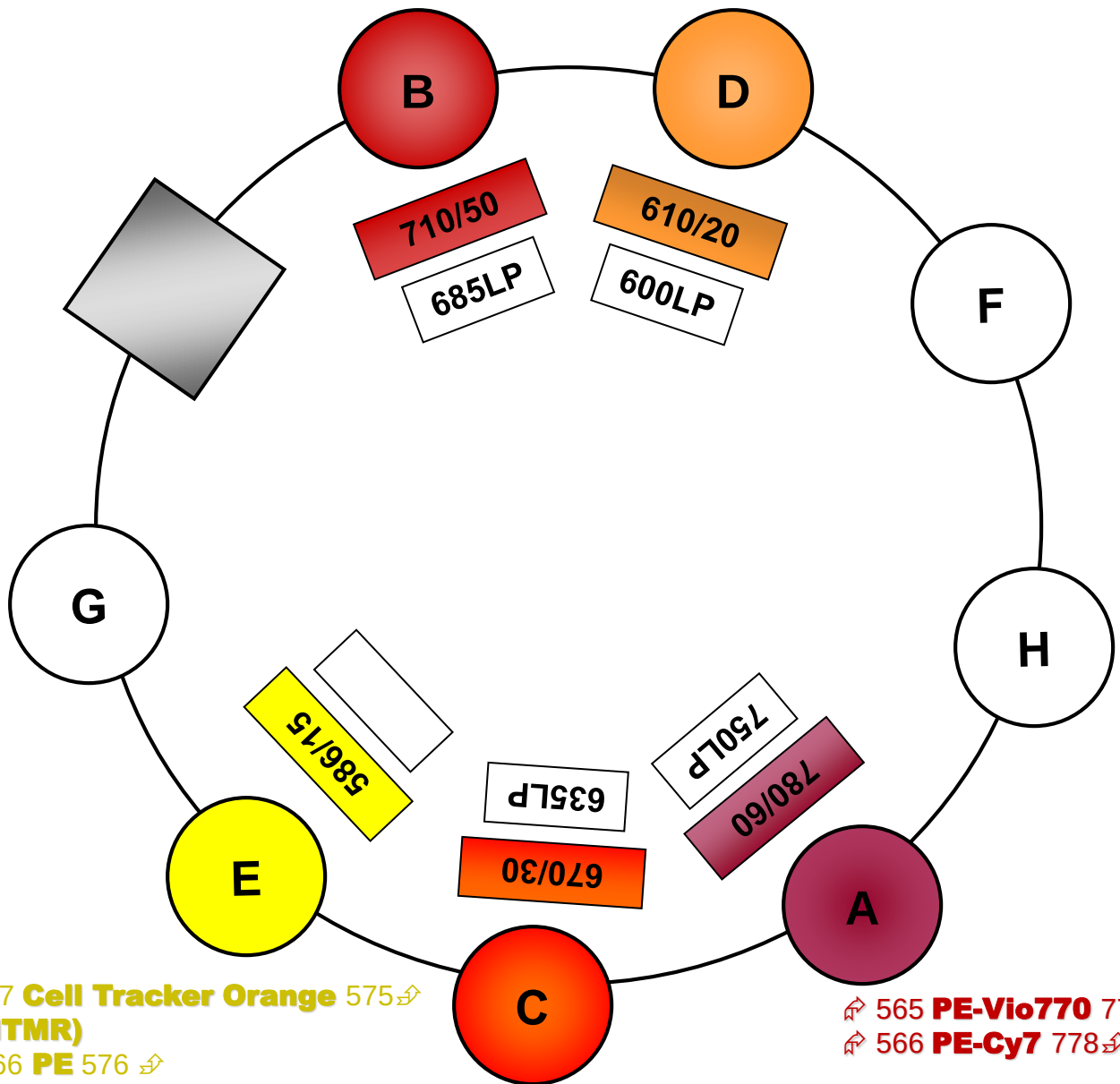

↻ 547 **Cell Tracker Orange (CMTMR)** 575 ↻

↻ 566 **PE** 576 ↻

↻ 553 **Ds-Red** 586 ↻

↻ 351 **PI** 617 ↻

↻ 565 **PE-Vio770** 775 ↻

↻ 566 **PE-Cy7** 778 ↻

↻ 543 **7-AAD** 647 ↻

↻ 566 **PE-Cy5** (↻ 647) 666 ↻

Configuration FORTESSA X20

Violet Laser (405 nm)

↻ 300 **Qdot 705** 702 ↻

↻ 405 **BV711** 711 ↻

↻ 414 **Super Bright 702** 702 ↻

↻ 300 **Qdot 605** 603 ↻ ↻ 405

BV605 603 ↻

↻ 405 **eVolve605** 610 ↻

↻ 414 **Super Bright 620** 620 ↻

↻ 405 **BV421** 421 ↻

↻ 400 **Zombie Violet FV423** ↻

↻ 398 **Cascade Blue** 420 ↻

↻ 401 **Alexa Fluor 405** 422 ↻

↻ 404 **Pacific Blue** 456 ↻

↻ 360 **Dapi Viab** 461 ↻


↻ 400 **VioBlue** 452 ↻

↻ 405 **eFluor450** 450 ↻

↻ 404 **BD Horizon V450** 448 ↻

↻ 406 **BD FVS 450** 450 ↻

↻ 414 **Super Bright 436** 436 ↻

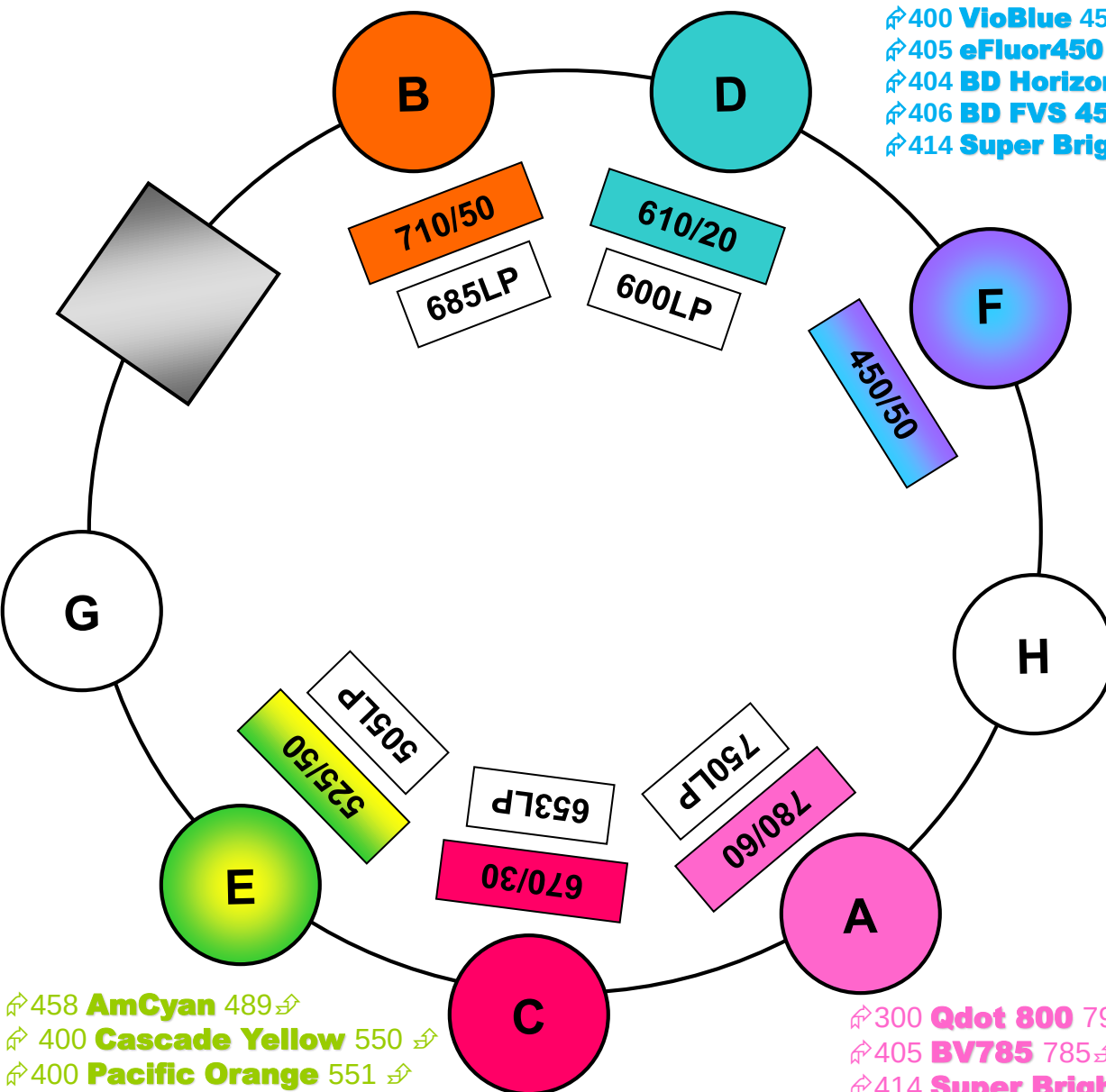

↻ 458 **AmCyan** 489 ↻

↻ 400 **Cascade Yellow** 550 ↻

↻ 400 **Pacific Orange** 551 ↻

↻ 405 **BV510** 510 ↻

↻ 415 **BD Horizon V500** 500 ↻

↻ 388 **VioGreen** 520 ↻

↻ 408 **BD FVS 510** 512 ↻

↻ 385 **Zombie Aqua FV** 515 ↻

↻ 300 **Qdot 800** 792 ↻

↻ 405 **BV785** 785 ↻

↻ 414 **Super Bright 780** 780 ↻

↻ 300 **Qdot 655** 654 ↻

↻ 405 **BV650** 645 ↻

↻ 405 **eVolve655** 655 ↻


↻ 414 **Super Bright 645** 645 ↻

Configuration FORTESSA X20

Red Laser (640 nm)

↪ 650 **APC-Alexa Fluor 750** (↪ 752) 779↪
↪ 650 **APC-H7** 780↪
↪ 650 **APC-Cy7** (↪ 753) 785↪
↪ 749 **Alexa Fluor 750** 775↪
↪ 652 **APC-Vio770** 775↪
↪ 650 **APC-eFluor780** 780↪
↪ 759 **BD FVS780** 780↪

↪ 650 **APC** 660↪
↪ 631 **Alexa Fluor 633** 647↪
↪ 630 **eFluor660** 660↪
↪ 650 **Alexa Fluor 647** 668↪
↪ 647 **Cy5** 660↪
↪ 649 **BD FVS 660** 660↪

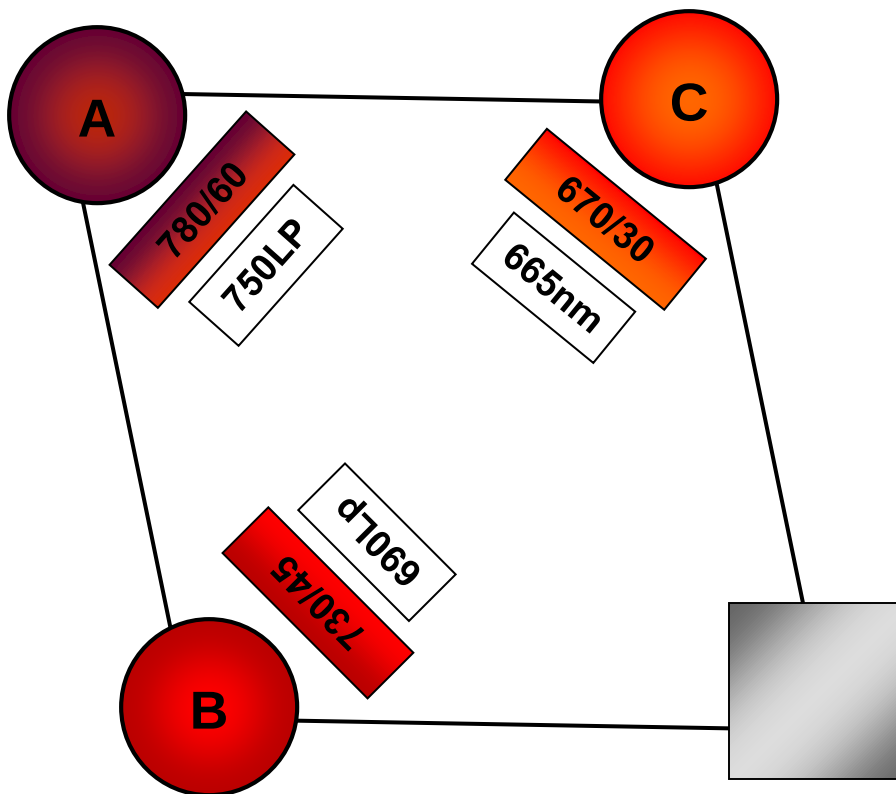

↪ 650 **APC-Cy5.5** (↪ 675) 695↪
↪ 675 **Cy5.5** 695↪
↪ 680 **Alexa Fluor 680** 700↪
↪ 700 **Alexa Fluor 700** 720↪
↪ 652 **APC-R700** 704↪
↪ 657 **BD FVS700** 700↪